


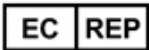













Doc. No. - REF.90.005

Revision: E

Effective Date: July 6, 2022

DCR: 22-100











SYMBOL	EXPLANATION OF REFERENCE SYMBOL	SYMBOL	EXPLANATION OF REFERENCE SYMBOL
	Catalog Number / "Reference" This symbol immediately precedes the part number of the product.		Batch Code / "Lot" This symbol immediately precedes the lot number of the product. G&H uses this for part traceability.
	Manufacturer Symbol for Manufacturer. May also include the date of manufacture in date format YYYY-MM-DD.		Authorized Representative in the European Community/European Union This symbol immediately precedes the name of our representative within the member states.
	Medical Device This symbol indicates the product is a medical device.		Manufacture Date Date of manufacture listed under or beside this symbol in YYYY-MM-DD format. <u>Note:</u> Not required if date of manufacture is used with Manufacturer symbol.
	CE Mark European Conformance to Medical Device Directive 93/42/EEC or European Conformance to Medical Device Regulation 2017/745	2797	Notified Body Identification Number BSI Group The Netherlands B.V. Say Building, John M. Keynesplein 9, 1066 EP Amsterdam, Netherlands
 (01) 10190886010433 (10) 123456 (11) 160929	UDI Datamatrix / "UDI Barcode" Manufacturer Unique Device Identifier "UDI" code used to contain traceability data of the part.		Unique Device Identifier Indicates the UDI carrier. <u>Note:</u> Used only when multiple UDI carriers are present on the label.
	R_x Only For use or distribution by orthodontic professionals only.		Do Not Reuse / "Single Use" Devices with this symbol are intended for single use.
	Chromium and/or Nickel warning Caution: Product contains Chromium and/or Nickel. To alert of the potential for allergic reactions to these known allergens.		Caution/Warning To bring awareness to the user of caution to prevent undesirable affects.
	Consult Instructions for Use / "IFU Book" For products with additional instructions inside the package.	 eIFU indicator	Consult electronic Instructions for Use For products with additional instructions electronically. <u>Note:</u> Web address will replace eIFU indicator text.

Doc. No. - REF.90.005

Revision: E

Effective Date: July 6, 2022

DCR: 22-100

SYMBOL	EXPLANATION OF REFERENCE SYMBOL	SYMBOL	EXPLANATION OF REFERENCE SYMBOL
	<p>Contains or presence of natural rubber latex Used on latex elastics.</p>		<p>Keep away from Sunlight/Heat Used on elastics. Product may demonstrate reduced performance / longevity when exposed to sunlight.</p>
	<p>Use-by-date The product should be used prior to the date located next to this symbol. Date format YYYY-MM-DD.</p>		<p>Temperature Limit Controlled storage is required. Number indicated to lower left is the min (lowest) temperature and number indicated to upper right is the max (highest) temperature at which the device should be stored in or exposed to.</p>
	<p>Country of Manufacture Identify the country of manufacture of the product. <u>Note:</u> CC will be replaced by the Alpha-3 country code.</p>		<p>Translation Indicates the original information has undergone translation. <u>Note:</u> Name and address of the entity that is responsible for translating activity shall be adjacent to the symbol.</p>
	<p>MR Safe Symbol used to identify product is safe to wear during MRI.</p>		<p>MR Conditional Symbol used to identify product is safe to wear during MRI but may affect imaging if performed above the neck.</p>
	<p>Hook Symbol used to identify brackets that have a hook.</p>		<p>Datamatrix / "Barcode" Manufacturer code used to contain traceability data of the part.</p>